BLACKROLL®


HEAT HEAD

Instruction Manual

EN | DE | FR | NL | IT | PL | CN

Product Overview

The BLACKROLL® HEAT HEAD is a heated massage attachment specially developed for the BLACKROLL® FASCIA GUN. It combines pleasant warmth with vibration – ideal for supporting relaxation and recovery before or after training.

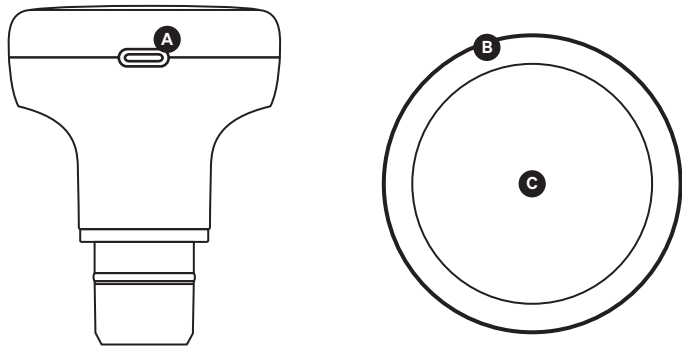Instructions for use

- A** USB-C port
- B** Indicator Light
- C** Metal Cover

Technical specifications

Compatibility:	Only BLACKROLL® FASCIA GUN
Operating mode:	Automatic activation via vibration
Trigger frequency:	from approx. 20 Hz (Level 1)
Battery:	600 mAh Li-Ion (2.2 Wh)
Runtime:	approx. 90 minutes
Charging port:	USB-C (Input 5V-1A)
Heating temperature:	max. 55 °C
Perceived temperature:	~45 °C
Weight:	approx. 63 g

First Use & Charging

Important: Before first use, connect the HEAT HEAD to a USB-C cable and charger to deactivate transport mode and fully charge the battery.

Charging

1. Connect the HEAT HEAD to a USB-C cable (not included).
2. The LED will blink while charging.
3. A solid LED indicates full charge.
4. Charging time: approx. 2 hours.

Note: Always fully charge before initial use.

Operation

1. Firmly attach the HEAT HEAD to your BLACKROLL® FASCIA GUN.
2. Switch on the gun – the head detects vibration automatically.
3. The heating function activates from approx. 20 Hz (Level 1).
4. Target temperature is reached after about 10–20 seconds.

Note: The heating function only remains active while the gun is vibrating.

Cleaning & Care

- + Clean with a soft, slightly damp cloth after each use
- + Do not use harsh cleaning agents
- + Do not immerse in water
- + Store in a cool, dry, and protected place

Safety Instructions / Contraindications

Do not use if:

- + You have open wounds, infections, or skin diseases in the treatment area
- + You have active tumors, fever, or acute inflammation
- + You have recent fractures, surgeries, or acute injuries
- + You wear pacemakers or other electronic implants
- + You suffer from thrombosis, blood clotting disorders, or high tendency to bruise
- + You are pregnant (unless medically approved)
- + Keep out of reach of children

Use only after medical consultation if:

- + You have rheumatism, fibromyalgia, or osteoporosis
- + You suffer from spinal disc or neurological conditions
- + You take anticoagulant medication
- + You have metal implants in the treatment area

Disposal & Environment

- + Clean with a soft, slightly damp cloth after each use
- + Do not use harsh cleaning agents
- + Do not immerse in water
- + Store in a cool, dry, and protected place

Warranty

The seller grants the buyer, provided the buyer is a consumer, the statutory warranty in accordance with legal regulations for a period of 2 years from the date of purchase of the item. The seller's warranty obligations are initially limited to supplementary performance, whereby the seller is free to choose whether to deliver a defect-free item or to repair the defective item. If the supplementary performance fails, the buyer retains the right to assert further statutory claims. The seller excludes warranty claims for damages that arise from accidents, obviously incorrect use, misuse, or improper transport; from repairs carried out by unauthorized personnel; from use not in accordance with the operating instructions; or from improper handling of the purchased item. Furthermore, claims for damages by the buyer are excluded. Exempt from this are the buyer's claims for damages resulting from injury to life, body, or health, as well as liability for other damages based on an intentional or grossly negligent breach of duty by the seller, its legal representatives, or agents. Liability also remains unaffected for breaches of essential contractual obligations whose fulfillment makes the proper execution of the contract possible in the first place and on the compliance with which the buyer may regularly rely. In the event of a slightly negligent breach of such contractual obligations, the seller is liable only for the foreseeable, contract-typical damage, unless it concerns claims for damages by the customer resulting from injury to life, body, or health.

Note: This product is not a medical device and does not replace medical treatment.

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BLACKROLL AG
WSVK Oederan GmbH